


**WARRANTY:** B. BRAUN INTERVENTIONAL SYSTEMS INC. WARRANTS THAT THIS PRODUCT WAS MANUFACTURED ACCORDING TO APPLICABLE STANDARDS AND SPECIFICATIONS. PATIENT CONDITION, CLINICAL TREATMENT, AND PRODUCT MAINTENANCE MAY AFFECT THE PERFORMANCE OF THIS PRODUCT. USE OF THIS PRODUCT SHOULD BE IN ACCORDANCE WITH THE INSTRUCTIONS PROVIDED AND AS DIRECTED BY THE PRESCRIBING PHYSICIAN.

# ACCEL Centesis Catheter

## Percutaneous Fluid Drainage Catheters

### INSTRUCTIONS FOR USE

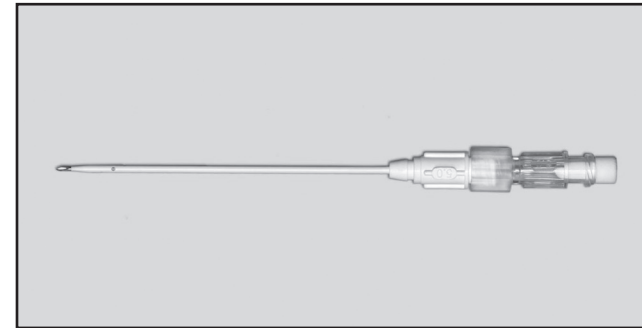

LS-00174-01-AA 2016-04

Rx only

Store at room temperature



Do not reuse



Do not resterilize



Date of manufacture



Not made with DEHP



Not made with natural rubber latex



Catalog number



Do not use if package is damaged



Sterilized using ethylene oxide



Keep away from sunlight



Use by



Batch code

*Interventional Systems*

**B | BRAUN**

Distributed by:

**B. Braun Interventional Systems Inc.**

824 Twelfth Avenue  
Bethlehem, PA 18018  
www.bisusa.org

Customer Service, ordering

TEL: (877) 836-2228

FAX: (610) 849-1334

Technical Support

TEL: (800) 443-VENA (8362)

Made in U.S.A.

*Interventional Systems*  
**B | BRAUN**

Distributed by:

**B. Braun Interventional Systems Inc.**

824 Twelfth Avenue  
Bethlehem, PA 18018  
www.bisusa.org



Manufacturer:  
pfm medical inc.  
1815 Aston Ave.,  
Suite 106  
Carlsbad, CA 92008

**INDICATIONS FOR USE:**

The ACCEL Centesis Catheter is for short term percutaneous fluid drainage.

**CONTRAINDICATIONS:**

This product should only be used by qualified personnel who are familiar with the technique. During insertion avoid contact with bone, cartilage, and scar tissue which can damage the catheter tip.

**CAUTION:**

- For Single Patient Use Only. Do not attempt to clean or resterilize this product.
- After use, this product may be a potential biohazard.
- Dispose in accordance with applicable laws and regulations.
- Handle in a manner that will prevent contamination.
- Do not use a device that has been damaged.
- Do not use any device if the package is opened or damaged.

**WARNING:**

This product contains chemicals known to the state of California to cause cancer and birth defects or other reproductive harm.

**NOTE:**

These instructions for the ACCEL Centesis Catheter are not meant to define or suggest any medical or surgical technique. The individual physician is responsible for the proper procedure and techniques to be used with this product.

French Size	Needle O.D.	Sheath Length	Catheter O.D.	Catheter I.D.
4	20g (.09 mm)	7cm	1.4mm	1.0mm
4	20g (.09 mm)	10cm	1.4mm	1.0mm
5	19g (1.1 mm)	7cm	1.7mm	1.2mm
5	19g (1.1 mm)	10cm	1.7mm	1.2mm
5	19g (1.1 mm)	15cm	1.7mm	1.2mm

**STORAGE:**

Store in a cool, dark, dry place.

**PROCEDURE**

( Single step catheterization through a skin incision )

1. Perform a skin incision under local anesthesia.
2. Introduce the ACCEL Centesis Catheter into the cavity using ultrasound guidance, CT or fluoroscopy.
3. Unlock the needle locking collar and remove the needle from the catheter.
4. Connect a syringe or suction tube to the catheter hub.

**ORDERING INFORMATION**

Reference Number	Description	Quantity
613210	ACCEL Centesis Catheter 4Fr, 7cm Length	5 per pack
613211	ACCEL Centesis Catheter 4Fr, 10cm Length	5 per pack
613212	ACCEL Centesis Catheter 5Fr, 7cm Length	5 per pack
613213	ACCEL Centesis Catheter 5Fr, 10cm Length	5 per pack
613214	ACCEL Centesis Catheter 5Fr, 15cm Length	5 per pack