

PTS® Percutaneous Transluminal Sizing

Balloon Catheters



Catheter Characteristics

- The PTS® Sizing Balloon catheter has been designed for use in patients with cardiovascular defects where accurate measurement of the defect is important to select the appropriately sized occluder device.

Super Thin Balloon

- The PTS® Sizing Balloon catheter is super-thin for a low deflated profile that maintains tip flexibility.


Double Tapered Balloon

- The PTS® Sizing Balloon catheter has short tapers at the distal and proximal ends for enhanced transition across the stenosis and post dilatation ease of withdrawal into the introducer.


Designed for controlled low pressure inflation

Tip Details


3.0 cm Balloon Length

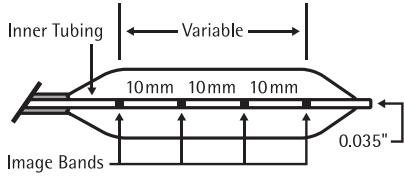

4.0 cm Balloon Length

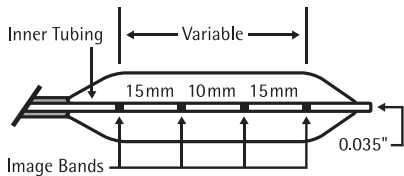

5.0 cm Balloon Length

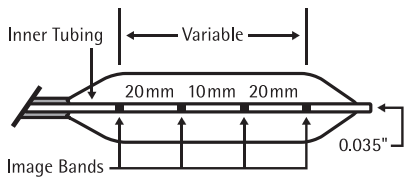

Ordering Information:

For more information or to place an order, contact your B. Braun Interventional Systems representative or call 1-877-VENA CAV (836-2228)

Item No.	Balloon Diameter (mm)	Balloon Length (cm)	Introducer Size (FR)	Shaft Size (FR)	Usable Length (cm)	Guide Wire (inches)	Rated Burst (ATM)	Catalog No.
613700	20.0	3.0	8	8	80	0.035	1.5	PTS203
613701	25.0	3.0	8	8	80	0.035	1.5	PTS253
613702	30.0	3.0	8	8	80	0.035	1.0	PTS303
613703	30.0	4.0	8	8	80	0.035	1.0	PTS304
613704	30.0	5.0	8	8	80	0.035	1.0	PTS305
613705	40.0	3.0	9	8	80	0.035	0.5	PTS403
613706	40.0	4.0	9	8	80	0.035	0.5	PTS404
613707	40.0	5.0	9	8	80	0.035	0.5	PTS405

The PTS® Sizing Balloon Catheters are recommended for use in those patients with cardiovascular defects wherein accurate measurement of the defect is important to select the appropriately sized occluder device.

Refer to the PTS® Sizing Balloon instructions for use for relevant warnings, precautions and complications.

Caution: Federal law restricts this device to sale by or on the order of a physician.

Rx only

CV1050-1/11
PTS® is a registered trademark of NuMED, Inc.
©2011 B. Braun Interventional Systems Inc.



Components and packaging do not contain natural rubber latex



Manufactured by:
NuMED, Inc.
2880 Main St.
Hopkinton, NY 12965
315-328-4491



Manufactured for:
B. Braun Interventional Systems Inc.
824 Twelfth Avenue
Bethlehem, PA 18018
Tel: 1-877-VENA CAV (836-2228)
Fax: 610-266-3982
www.bbabraunusa.com